






BD Biliary

ELLA-BD Stent® Biliary THP

**Worldwide unique self-expandable
biodegradable biliary stent**

For management of benign
biliary strictures

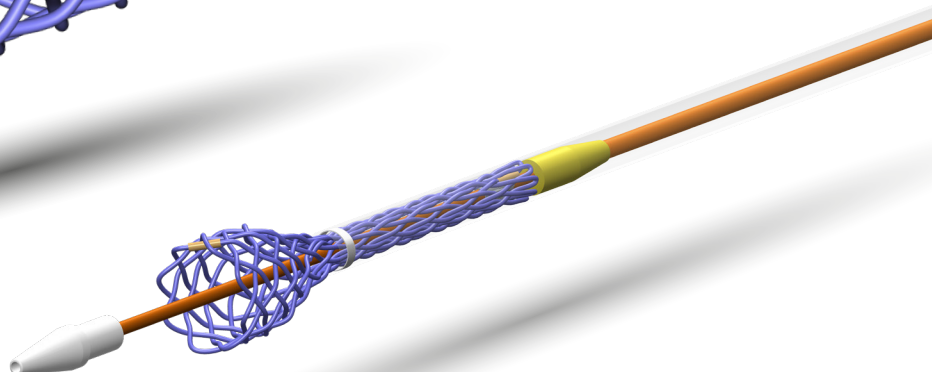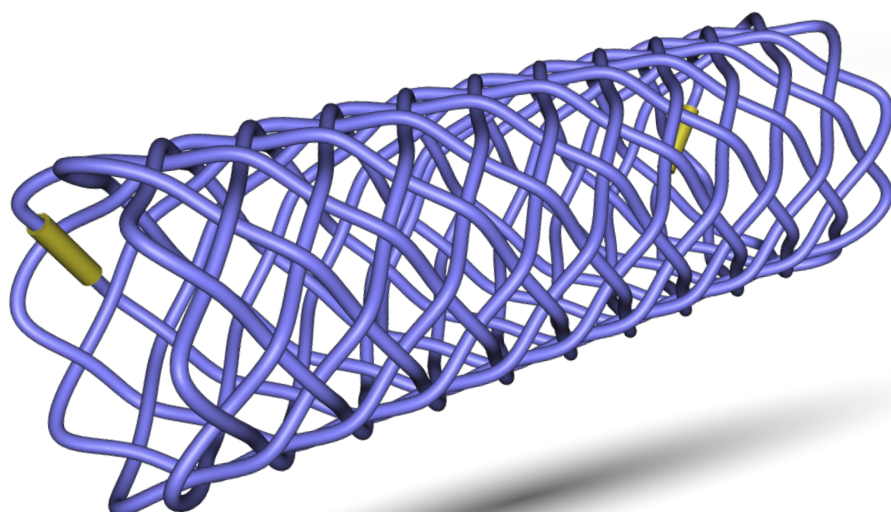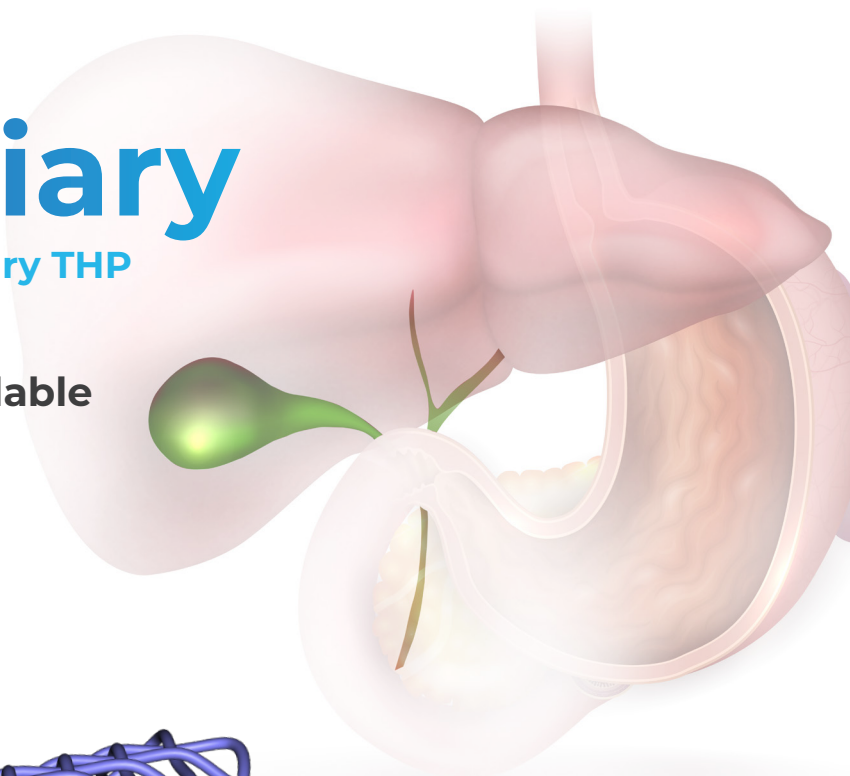

Designed & manufactured
in the Czech Republic (EU)

3D ANIMATION



www.ellacs.eu

Manufacturer:
ELLA-CS, s.r.o.
Milady Horakove 504/45, Trebes
500 06 Hradec Kralove, Czech Republic

Indication

ELLA-BD Stent® Biliary THP is indicated for the treatment of adult patients (≥ 18 yrs) with benign biliary strictures.

Technical information

Delivery system 11.8 F available in standard length of 60 cm. Simple stent compression procedure – compression tool included.



Recommended guide wire: 0.035" (0.89 mm) ultra stiff / Minimum length 180 cm

Available sizes

REF No.	Nominal stent body diameter [mm]	Nominal stent length [mm]	Delivery system
019-15-08-030	8	30	11.8 F 60 cm
019-15-08-040		40	
019-15-08-050		50	
019-15-08-060		60	
019-15-08-070		70	
019-15-08-080		80	
019-15-10-030	10	30	
019-15-10-040		40	
019-15-10-050		50	
019-15-10-060		60	
019-15-10-070		70	
019-15-10-080		80	

Features and Benefits

- + **Biodegradable material** – degrades usually within 3-6 months after the implantation, no need for removal – less hospital visits, lower clinical risks.
- + **Self-expandable** – stent integrity and radial force are maintained for 6 to 8 weeks after implantation.
- + **Uncovered design** – low migration.
- + **Radiopaque markers** – accurate stent positioning.

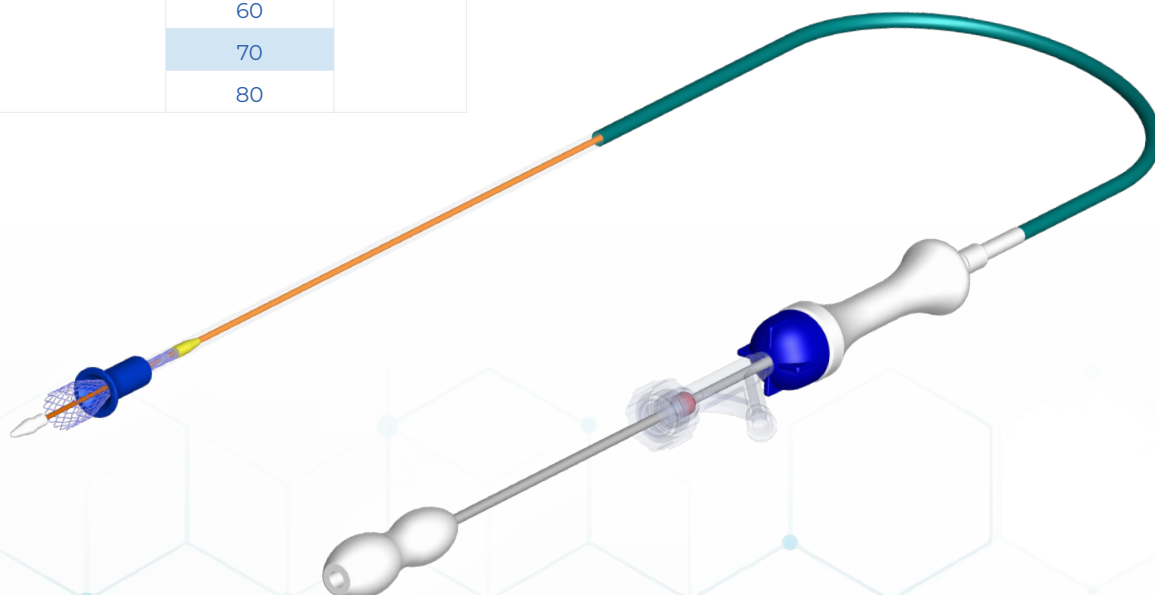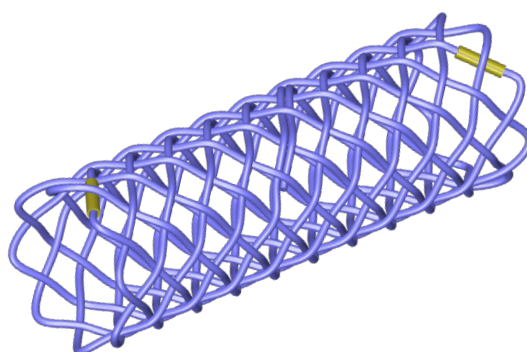

PUBLICATIONS



www.ellacs.eu

info@ellacs.eu

CE 2797

Manufacturer:
ELLA-CS, s.r.o.
Milady Horakove 504/45, Trebes
500 06 Hradec Kralove, Czech Republic

SX-15-AJ-02/24/P/REV-0-02/2024